

FOREO

LUNA™ 3 plus

**THERMO-FACIAL
CLEANSING
& FIRMING MASSAGE**

plus Microcurrent



FULL USER MANUAL



GETTING STARTED

Congratulations on taking the first step toward discovering smarter skincare by acquiring LUNA 3 plus. Before you begin to enjoy all the benefits of sophisticated skincare technology in the comfort of your home, please take a few moments to carefully read the instructions in this manual.

Please READ ALL INSTRUCTIONS BEFORE USE and utilize this product only for its intended use as described in this manual.

⚠ WARNING: NO MODIFICATION OF THIS EQUIPMENT IS ALLOWED.

LUNA 3 plus OVERVIEW

Put the power of beauty into your own hands with the only beauty-tech in the world that features heated cleansing and targeted microcurrent in one smart device. Featuring an indulgent skincare ritual, this ultra-soft heated cleansing device melts away impurities while targeted microcurrent restores the skin's youthful radiance. With LUNA 3 plus clean, younger-looking skin has never been more attainable.

GETTING TO KNOW LUNA 3 plus

T-SONIC PULSATIONS

Effectively cleanse with up to 8,000 pulsations per minute at 16 different speeds

THERMO-TECH

Thermal touchpoints warm skin to the perfect temperature for a deep yet gentle cleanse

EXTRA-SOFT & HYGIENIC SILICONE

35x more hygienic than nylon bristles and suitable for all skin types. BPA & phthalate free

MICROCURRENT PINS

Deliver microcurrent directly to skin for non-invasive age renewal

MICROCURRENT MODE

Targets signs of aging for firm, younger-looking skin

UNIVERSAL POWER BUTTON

Turns device on and off and puts device in Bluetooth pairing mode

INDICATOR LIGHT

Flashes to indicate when device is in Bluetooth pairing mode or when device needs to be charged

CHARGING PORT

Up to 125 uses per charge. 100% waterproof closure



2ML SAMPLE SERUM

Use with microcurrent mode for optimal conductivity and enhanced results



INTERACTIVE FOREO APP

Control device settings, cleanse intensity & access microcurrent treatments



USB CABLE

Charge anytime, anywhere



TRAVEL POUCH

Conveniently stores device for skincare anywhere

HOW TO USE LUNA 3 plus

CLEANSING MODE

To unlock and register your device for first time use, activate Bluetooth on your smartphone and press the universal button on your LUNA 3 plus to sync to the FOREO app. Next, select your cleansing preferences. We recommend cleansing with LUNA 3 plus daily for a refined, healthy glow.

- 1 Remove all makeup, dampen skin and then apply the FOREO Micro-Foam Cleanser.
- 2 Press universal button again to activate cleansing mode.
- 3 Cleanse using circular movements on cheeks and forehead, gliding up and down the nose.
- 4 Press and hold button for 3 seconds to turn off device.
- 5 Rinse and pat your face dry. Apply other skincare products, if any.



MICROCURRENT MODE

- 1 Start with clean, dry skin.
- 2 Apply FOREO SERUM SERUM SERUM to treatment area.
- 3 Follow instructions in the FOREO app to activate your preferred microcurrent treatment.
- 4 Follow with moisturizer of your choice.

** Device will automatically switch off after routine is complete.*

CAUTION: LUNA 3 plus must be completely dry before activating your microcurrent treatment for safe and effective use.



CLEANING YOUR DEVICE

Always clean your device thoroughly after use. Wash the brush surface with water and soap, then rinse with warm water. Avoid using clay-based, silicone-based or grainy cleansers, exfoliators or scrubs, as they may damage the device's soft silicone touchpoints. Pat dry with a lint-free cloth or towel. After use, we recommend spraying the device with FOREO Silicone Cleaning Spray and rinsing with warm water for optimal results.

NOTE: Never use cleaning products containing alcohol, petrol or acetone, as they may irritate the skin and damage the silicone.

IMPORTANT FOR OPTIMUM SAFETY:

- If you have a skin condition or any medical concerns, please consult a physician before use.
- Do not use if you have any existing medical conditions such as epilepsy, a hemorrhagic disease, cancer, tumors, or perceptual disturbance.
- Do not use if you have any implanted medical device or any other electronic instrument or body aids.
- Do not use if you have had a laser treatment, chemical peel, or have any other form of injured or damaged skin.
- Do not use the microcurrent feature on the chest / breast area, eye area (circular muscle within the orbital rim), mid-line (bone) of neck, or the genitals / groin area.
- Do not use if you have had cosmetic surgery on your face.
- If you have suspected or diagnosed heart disease, follow precautions recommended by your physician.
- The long-term effects of microcurrent are unknown.
- The safety of the device for use during pregnancy has not been established. If you are pregnant, consult your physician before using this product.
- There may be a perception of flashing lights during the use of your device, resulting from the stimulation of your optic nerve. If you notice this condition constantly without the use of your device, consult your physician.
- Electronic monitoring equipment, such as ECG monitors and ECG alarms, may not operate properly when the device is in use.
- Cleansing with LUNA 3 plus should be comfortable – if you experience discomfort, discontinue use immediately and consult a physician.
- Exercise particular care when cleansing the under-eye regions and do not bring the device into contact with the eyelids or eyes themselves.
- For reasons of hygiene, we do not recommend sharing your LUNA 3 plus with anyone else.
- Avoid leaving your LUNA 3 plus in direct sunlight and never expose it to extreme heat or boiling water.
- Close supervision is necessary when this device is used by, on, or near children, as well as those with reduced physical and mental abilities.
- Discontinue use if this product appears damaged in any way. This product contains no serviceable parts.
- Given the efficiency of the FOREO cleansing routine, we recommend that you do not use LUNA 3 plus for more than 3 minutes at a time when in cleansing mode.
- Do not use the device in microcurrent mode if it has been submerged in water. LUNA 3 plus should be completely dry before activating the microcurrent function.
- A slight tingling sensation may occur while using the device. Lowering the intensity may reduce or eliminate this sensation.
- The device is intended for cosmetic use on the face. Any harmful consequences resulting from misuse, application to other body areas, connection to improper voltage sources, dirty conductive solution or spheres, or any other improper applications are not the responsibility of FOREO.
- Use this device only for its intended use as described in this manual. If you do not find the answer to your specific question, or if you have any other questions regarding the device's operation, please visit www.foreo.com

TROUBLESHOOTING

PRECAUTIONS TO BE TAKEN IN THE EVENT OF CHANGES IN THE PERFORMANCE OF LUNA 3 PLUS.

If LUNA 3 plus is not activated when pressing the center button:

- Battery is empty. Recharge using USB charging cable for up to 4 hours until fully charged and then restart your device by holding down the power button.

If LUNA 3 plus cannot be switched off and/or interface buttons do not respond:

- Microprocessor is temporarily malfunctioning. Press and hold power button to restart the device.

If LUNA 3 plus won't sync to the FOREO app:

- Switch your Bluetooth off and then on again to try reconnecting.
- Close the FOREO app and then reopen it to start the process over.
- Check to see if the app needs to be updated in your mobile device's app store.

WARRANTY TERMS & CONDITIONS

REGISTER WARRANTY

To activate your 2-Year Limited Warranty register through the FOREO app, or visit foreo.com/product-registration for more information.

2-YEAR LIMITED WARRANTY

FOREO warrants this device for a period of TWO (2) YEARS after the original date of purchase against defects due to faulty workmanship or materials arising from Normal Use of the device. The warranty covers working parts that affect the function of the device. It does NOT cover cosmetic deterioration caused by fair wear and tear, or damage caused by accident, misuse or neglect. Any attempt to open or take apart the device (or its accessories) will void the warranty.

If you discover a defect and notify FOREO during the warranty period, FOREO will, at its discretion, replace the device free of charge. Claims under warranty must be supported by reasonable evidence that the date of the claim is within the warranty period. To validate your warranty, please keep your original purchase receipt together with these warranty conditions for the duration of the warranty period.

To claim your warranty, you must log in to your account at www.foreo.com and then select the option to make a warranty claim. Shipping costs are nonrefundable. This undertaking is in addition to your statutory rights as a consumer and does not affect those rights in any way.

DISPOSAL INFORMATION

Disposal of old electronic equipment (applicable in the EU and other European countries with separate waste collection systems).




The crossed-out dustbin symbol indicates that this device should not be treated as household waste, but rather be brought to the appropriate collection point for recycling of electrical and electronic equipment. By ensuring this device is disposed of correctly, you will help prevent potentially negative consequences for the environment and human health which could be caused by inappropriate waste handling of the product. The recycling of materials will also help conserve natural resources.

For more detailed information about the recycling of your device, please contact your local household waste disposal service or your place of purchase.

BATTERY REMOVAL

NOTE: This process is not reversible. Opening the device will void its warranty. This action must only be undertaken when the device is ready to be disposed of.

Because this device contains a lithium-ion battery, the battery must be removed before disposal and should not be thrown away with household waste. To remove the battery, open the inner plastic shell after removing silicone outer layer and remove the battery to be disposed of in accordance with your local environmental regulations. Wear gloves during this process for your safety. Detailed visual instructions are provided below:

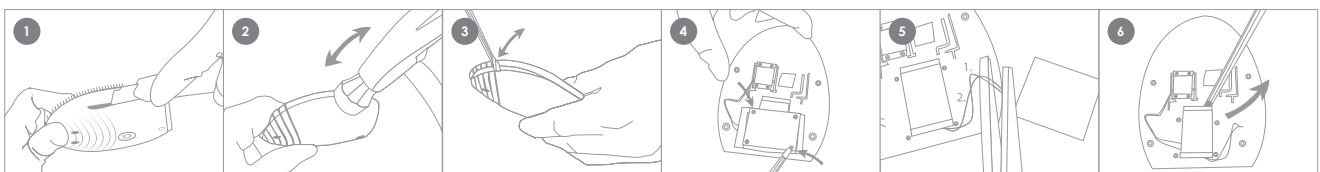

SPECIFICATIONS

MATERIALS: Body-safe silicone, ABS plastic, polycarbonate, zinc alloy, aluminium & gold

COLOR: Pink & Lavender

SIZE: 102*39.2 *82.5mm

WEIGHT: 155g

BATTERY: Li-ion 2000mAh 3.7V

USAGE: Up to 125

STANDBY: 90 days

FREQUENCY: 100 Hz

MAX NOISE LEVEL: <50dB

INTERFACE: 1- button

MICROCURRENT LEVELS: 15

SPEED LEVEL: 16

Disclaimer: Users of this device do so at their own risk. Neither FOREO nor its retailers assume any responsibility or liability for any injuries or damages, physical or otherwise, resulting, directly or indirectly, from the use of this device. Further, FOREO reserves the right to revise this publication and to make changes from time to time in the contents thereof without obligation to notify any person of such revision or changes.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

RF exposure statement: This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

©2021 FOREO AB. ALL RIGHTS RESERVED.

EU DISTR.: FOREO AB, KARLAVÄGEN 41, 114 31 STOCKHOLM, SWEDEN.

US DISTR.: FOREO INC., 1525 E PAMA LANE, LAS VEGAS, NV 89119, USA.

AUS DISTR.: 6 PROHASKY STREET, PORT MELBOURNE, VIC, 3207 AUSTRALIA.

DESIGNED AND DEVELOPED BY FOREO SWEDEN.

MANUFACTURED FOR FOREO AB.

WWW.FOREO.COM