

FOREO
FOR EVERYONE



FLIP™ play
advanced LED light hairbrush



ONLINE MANUAL

FULL USER MANUAL

GETTING STARTED

Congratulations on discovering smarter haircare with FLIP™ play advanced LED light hairbrush. Before you begin to enjoy all the benefits of sophisticated technology in the comfort of your home, please take a few moments to carefully read the instructions in this manual.

Please READ ALL INSTRUCTIONS BEFORE USE and utilize this product only for its intended use as described in this manual.

⚠ **WARNING:** NO MODIFICATION OF THIS EQUIPMENT IS ALLOWED.

FLIP™ play advanced LED hairbrush OVERVIEW

Introducing the hairbrush designed to improve hair & scalp health. With LED light, this innovative hairbrush detangles like a regular brush, while also supporting your natural hair cycle to help reduce breakage & thinning.

Red light helps revitalize follicles to boost thicker, healthier-looking hair. Blue light combats oiliness & dead skin to support a balanced scalp. Green light supports scalp & hair health to maintain healthy hair follicles.

Designed for all hair types. Yes - we do mean ALL.

GETTING TO KNOW FLIP™ play advanced LED hairbrush

OUTER FRAME

Flips forward to shield bristles from picking up lint / getting tangled in your bag. Flips back to open up bristles, ready to brush and repair your hair.

LED LIGHT THERAPY

Strengthens hair, balances scalp and supports natural hair cycle to combat hair loss.

ULTRA-FLEXIBLE BRISTLES

Tough on tangles, gentle on hair and scalp.

INDICATOR LIGHT


Flashes to indicate when device is in Bluetooth pairing mode, Smart Control mode, and when device needs to be charged.

UNIVERSAL POWER BUTTON

Turns device on and off, switches between LED colors, and puts device in Bluetooth pairing mode.

USB CHARGING PORT

100% waterproof. Up to 6 hours of use per charge, with included USB charging cable.

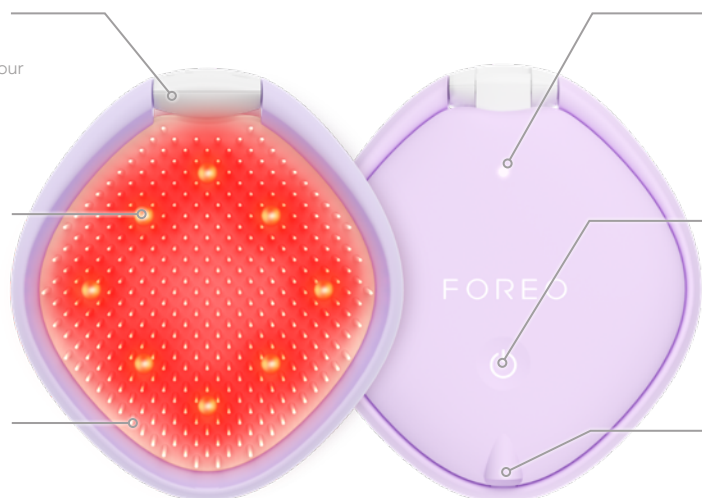

HOW TO USE FLIP™ play advanced LED hairbrush

UNLOCK YOUR DEVICE

Before using your hairbrush for the first time, you must unlock it. You can do this automatically by registering it via the FOREO app.

1. Download the FOREO For You app on your mobile phone by scanning the QR code.
2. Log in to your existing account, or sign up for a new one.
3. Add device (on the top of your screen).
4. Choose the device series.
5. Follow the instructions on the app and fill in the purchase information.

And your hairbrush is ready to use!



HOW TO USE



STEP 1

Flip open the brush, and press the universal button to turn on the LED light.



STEP 2

Press the universal button again to switch between LED light colors.



STEP 3

Brush hair as normal, from scalp to ends of hair. Once finished, press the universal button for 3 seconds to turn the LED light off.

IMPORTANT

FOR OPTIMUM SAFETY:

- Brushing your hair with FLIP™ play advanced LED light hairbrush should be comfortable – if you experience any discomfort, discontinue use immediately and consult a physician.
- If taking a photosensitizing medication or herbal remedy, seek medical advice prior to application of the appliance as use may harm the skin.
- Do not stare directly at the LED lights.
- This appliance is not intended for the treatment of medical conditions. Medical advice should be sought to address, for example, moles, skin rash, itchy skin, skin fungus or infection, skin bumps, or skin tags; a list of possible side effects, including blistering, pain, and scarring.
- This product should not be used by children. Close supervision is necessary when this product is used by, on, or near children and those with reduced physical and mental abilities.
- For reasons of hygiene, we do not recommend sharing your device with anyone else.
- Avoid leaving your device in direct sunlight and never expose it to extreme heat or boiling water.
- Discontinue use if this product appears damaged in any way.
- This product contains no serviceable parts.
- Do not tamper or attempt to disassemble the product as this may cause damage to the unit.
- Do not use while driving or operating heavy machinery.
- Do not insert any object into any opening of the device.
- Do not use the device if it is overheating or you suspect that it is malfunctioning.
- This device should only be used with a SELV power adapter.
- Before charging, make sure that the plug and socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.
- Do not use the device while charging. Discontinue use if this device or charger is not working properly or appears damaged in any way. Only use the power cord supplied with your device.
- The device contains a non-replaceable battery. The battery must be removed from the appliance before it is disposed of. The appliance must be disconnected from the supply mains when removing the battery, and the battery is to be disposed of safely.
- Use this device only for its intended use as described in this manual. If you do not find the answer to your specific question, or if you have any other questions regarding the device's operation, please visit foreo.com/support.

CLEANING FLIP™ play advanced LED hairbrush

We recommend cleaning your device once a week. Wash the brush surface with water and soap, then rinse with warm water. Pat dry with a lint-free cloth or towel.

NOTE: Never use cleaning products containing alcohol, petrol, or acetone, as they may irritate the scalp and damage the silicone.

TROUBLESHOOTING

Precautions to be taken in the event of changes in the performance of your FLIP™ play advanced LED light hairbrush:

If the hairbrush cannot be switched on/off, and the universal button does not respond:

- Battery is empty. Recharge using the USB charging cable for up to 1.5 hours until fully charged, and then restart your device by pressing the universal button.
- Microprocessor is temporarily malfunctioning. Press and hold the universal button to restart the device.
- If this is your first time using the hairbrush, it's locked. Follow the UNLOCK YOUR DEVICE instructions under the How To Use section of this manual.

If the hairbrush won't sync to the FOREO app:

- Ensure that the hairbrush is fully charged, and that your Bluetooth is turned on.
- Switch your Bluetooth off and then on again to try reconnecting.
- Close the FOREO app and then reopen it to restart the process.
- The app may need to be updated. Delete and reinstall the app, then make sure to restart your phone.
- Delete your phone's cookies and cache memory.

DISPOSAL INFORMATION

Disposal of old electronic equipment (applicable in the EU and other European countries with separate waste collection systems).




This device should not be treated as household waste, but rather be brought to the appropriate collection point for recycling of electrical and electronic equipment. By ensuring this device is disposed of correctly, you will help prevent the potential negative consequences for the environment and human health which could be caused by inappropriate waste handling of the device. The recycling of materials will also help conserve natural resources.

For more information about the recycling of your device, please contact your local household waste disposal service or your place of purchase.

BATTERY REMOVAL

NOTE: This process is not reversible. This action must only be undertaken when the device is ready to be disposed of.

Because this device contains a lithium-ion battery, the battery must be removed before disposal and should not be thrown away with household waste. In order to do this, remove the silicone outer layer, open the inner plastic shell, remove the battery, and dispose of it in accordance with your local environmental regulations. Wear gloves during this process for your safety. Detailed visual instructions are provided below:

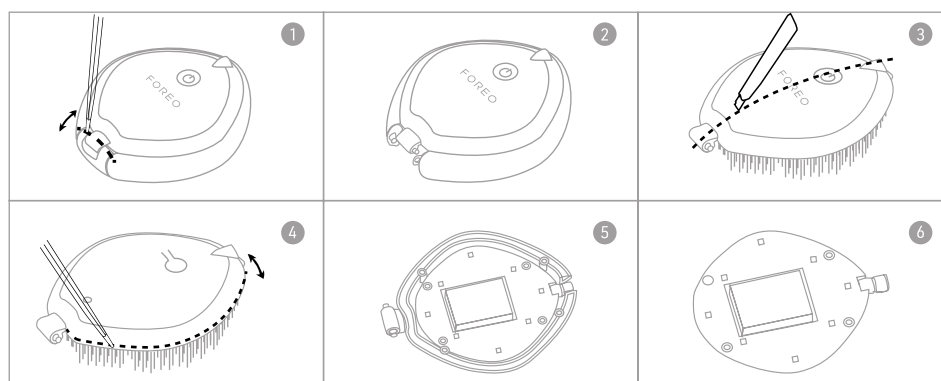

SPECIFICATIONS

MATERIALS: Body-safe silicone, PC + ABS plastic, TPEE

COLOR: I Lilac You! / Ice Ice Baby! / Midnight Magic!

SIZE: H101 x W92 x D26 mm

WEIGHT: 81.5 g

BATTERY: Li-ion 300mAh 3.7V

USES PER CHARGE: Up to 6 hours of use per 1.5 hour charge

STANDBY: 90 days

DISCLAIMER: Users of this device do so at their own risk. Neither FOREO nor its retailers assume any responsibility or liability for any injuries or damages, physical or otherwise, resulting directly or indirectly from the use of this device. Furthermore, FOREO reserves the right to revise this publication and to make changes from time to time in the contents thereof without obligation to notify any person of such revision or changes.

CAUTION: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This device complies with Part 15 of the FCC Rules and contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

(1) This device may not cause interference.

(2) This device must accept any interference, including interference that may cause undesired operation of the device.

FOREO declares that this device is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity for the product is available at www.foreo.com/support/eu-conformity.

Model may be changed for improvements without notice.

©2025 FOREO AB. ALL RIGHTS RESERVED. MANUFACTURED FOR FOREO AB.

EU DISTR.: FAQ EUROPE D.O.O., RADNICKA CESTA 202, 10 000, ZAGREB, CROATIA.

US DISTR.: FOREO INC., 1525 E PAMA LANE, LAS VEGAS, NV 89119, USA.

AUS DISTR.: FOREO AUS., 6 PROHASKY STREET, PORT MELBOURNE, VIC, 3207 AUSTRALIA.

UK RP: FOREO LTD., LAWFORD HOUSE, ALBERT PLACE, N3 1QA LONDON.

WWW.FOREO.COM